

Bed Panel Size	(2015x 540) ± 10mm
Bed height from the ground(Without pads)	Minimum height: (480 ± 20) mm Maximum height: (840 ± 20) mm
Lifting route (electric)	(340 ± 20) mm
Horizontal route (electric)	(350 ± 20) mm
Head and foot tilt (electric)	Head down tilt: 20° ± 2°
	Foot down tilt: 20° ± 2°
Lateral tilt (electric)	Left: 18° ± 2°
	Right: 18° ± 2°
Back plate fold (manual)	Upward: 70° ± 2°
	Downward: 32° ± 2°
Leg plate fold (manual)	Upward: 28° ± 2°
	Downward: 90° ± 2°
Head section fold (manual)	Upward: 30° ± 2°
	Downward: 90° ± 2°
Table panel X-ray absorption coefficient	≤ 0.44mmAL
Power supply	220V 50Hz 300VA
Net weight/ capacity weight	300kg / 400kg
Bed lock	Electric

## Standard Accessories



Y-001  
mattress pad




Y-002  
cable controller



Y-003  
lateral support



P-001/P-002  
lock block/clamping block

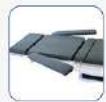

Y-004  
carbon fiber arm support



Y-006 fixed belt



Y-009 leg support



Y-008 anaesthesia screen



D-001 orthopedic tractor



D-004 hanging arm support



D-005 double deck arm support



D-006 exterior kidney bridge



D-007 knee crutches



D-002 cassette holder



D-008 adjustable spine pad



D-0010 upper limbs operation support



D-009 adjustable full carbon fibre spine shell

## DST-III Electric Hydraulic Integrated operating table



- 1.DST-III use electric and hydraulic. Germany Dunkerton motor with advanced hydraulic transmission.
- 2.350mm horizontal movement. Micro-computer control precision stability. Full carbon fiber material, China patent, very high clearance, supply accurate space for C-arm and X-ray imaging of full body without moving the patient.
- 3.304 stainless steel for the cover of ABS base. ≥100mm castors can be electric locked.
- 7.backup battery for one week (80-100times) surgeries.
- 8.It has standard and reverse standard for head and foot exchange function.
7. One reset function, with two independent electric controlling system. High density memory mattress.

Bed Panel Size	(2015x 540) ± 10mm
Bed height from the ground(Without pads)	Min height: (660 ± 20) mm Max height: (1000 ± 20) mm
Lifting route (electric)	(340 ± 20) mm
horizontal route(electric)	(350 ± 20) mm
Head and foot tilt (electric)	Head down tilt: 25° ± 2°
	Foot down tilt: 25° ± 2°
Lateral tilt (electric)	Left: 22° ± 2°
	Right: 22° ± 2°
Back plate fold(electric)	Upward: 70° ± 2°
	Downward: 32° ± 2°
Leg plate fold (manual)	Upward: 28° ± 2°
	Downward: 90° ± 2°
Head section fold (manual)	Upward: 30° ± 2°
	Downward: 90° ± 2°
Table panel X-ray absorption coefficient	≤ 0.44mmAL
Power supply	220V 50Hz 300VA
Net weight/ capacity weight	260kg / 300kg