



Knee Arthroscopy Surgery

Recommended accessories:

- 230000760 donut head gel pad x1
- 230000769 Torso Section gel pad x1
- 121000294 Arm rest x2
- 121000532 Knee support



Rectal Surgery

Recommended Accessories

- 121000294 Arm rest x2
- 121000347 Foot plate 1pair
- 121000348 Notched Clamp x2
- 121000355 Iliac support x1



Orthopedic Procedures

Recommended Accessories:

- 230000760 donut head gel pad x1
- 230000769 Torso Section gel pad x1
- 121000358 Orthopedics traction x1
- 121000349 High position arm holder x1
- 121000294 Arm rest x1



Thoracic Kidney/ Gall Bladder Surgery

Recommended Accessories

- 230000760 donut head gel pad x1
- 230000769 Torso Section gel pad x1
- 121000294 Arm rest x1
- 121000349 High position arm holder x1
- 230000766 Universal Dome Shaped Pad x1



Aegistab OP830/330

Electronic Operating Table



AEOMED
Reliable Quality Thoughtful Service

HQ: 11B2, Fengtai Science Park, Beijing, China (100070)
Science Park: No. 10, Chaobai Street, Yanjiao Development Zone, Sanhe City, Hebei Province, China (065201)
TEL: +86-10-5841 1198 FAX: +86-10-6371 8989
Http://www.aeonmed.com E-mail: marketing@aeonmed.com

The information contained in this leaflet is correct at the date of publication. The policy of Beijing Aeonmed Co., Ltd. is one of the continued improvement to its products. Because of this policy, Aeonmed reserves the right to make any change, which may affect the information in this leaflet without giving prior notices.
AOP830/330-1708

AEOMED
Reliable Quality Thoughtful Service

A cost-effective, multi-performance operating table holding up to 297Kg loading, adaptable to meeting different surgery requirement.



② High Efficiency Built-in Battery

- Imported from Panasonic
- Maximum 15 days working time
- No request for AC power to increase surgery safety



③ Perfect X-ray Access

- Remarkable 300mm longitudinal shifting capability create perfect X-ray access



④ Ergonomic Adjustment

- Convenient design for wrench of leg plate and head plate, move leg and head plate by one hand
- Imported gas springs guarantee smooth movement
- Imported round-hole clamps, easy to lock and unlock accessories, safe time during operation
- Adjustable table height (680-950mm) provides comfortable working condition



① Premium Quality Material

- Table frame is specially coated
- Other exposed steel is stainless steel SUS304
- Easy to clean and disinfect
- Optional with SUS304 Stainless steel table frame




⑤ Surgical Configuration

- Flex and re-flex function
- Table angles and section lengths can be adjusted within a wide range. 6 standard accessories and 13 optional ones
- Smart combination of all these meets various surgical needs.




⑥ Safety

- Control panel power off automatically when no operation over 1 minute, avoid any mis-operation

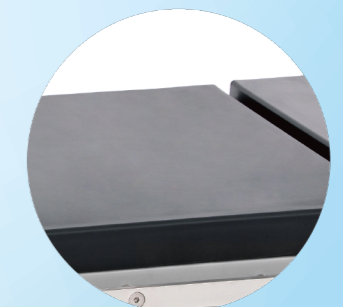


⑦ Pad

- PU Foam mattress, easy for disinfection and sterilization
- Effectively reduce possibility of patient bedsores

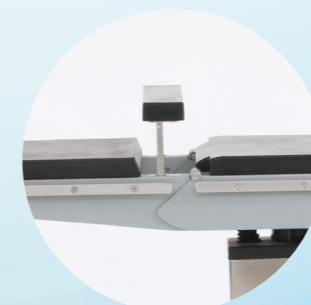


⑧ Built-in Elevator Kidney-bridge

- No negative effect on X-ray radiolucency, support kidney and gallbladder surgery efficiently


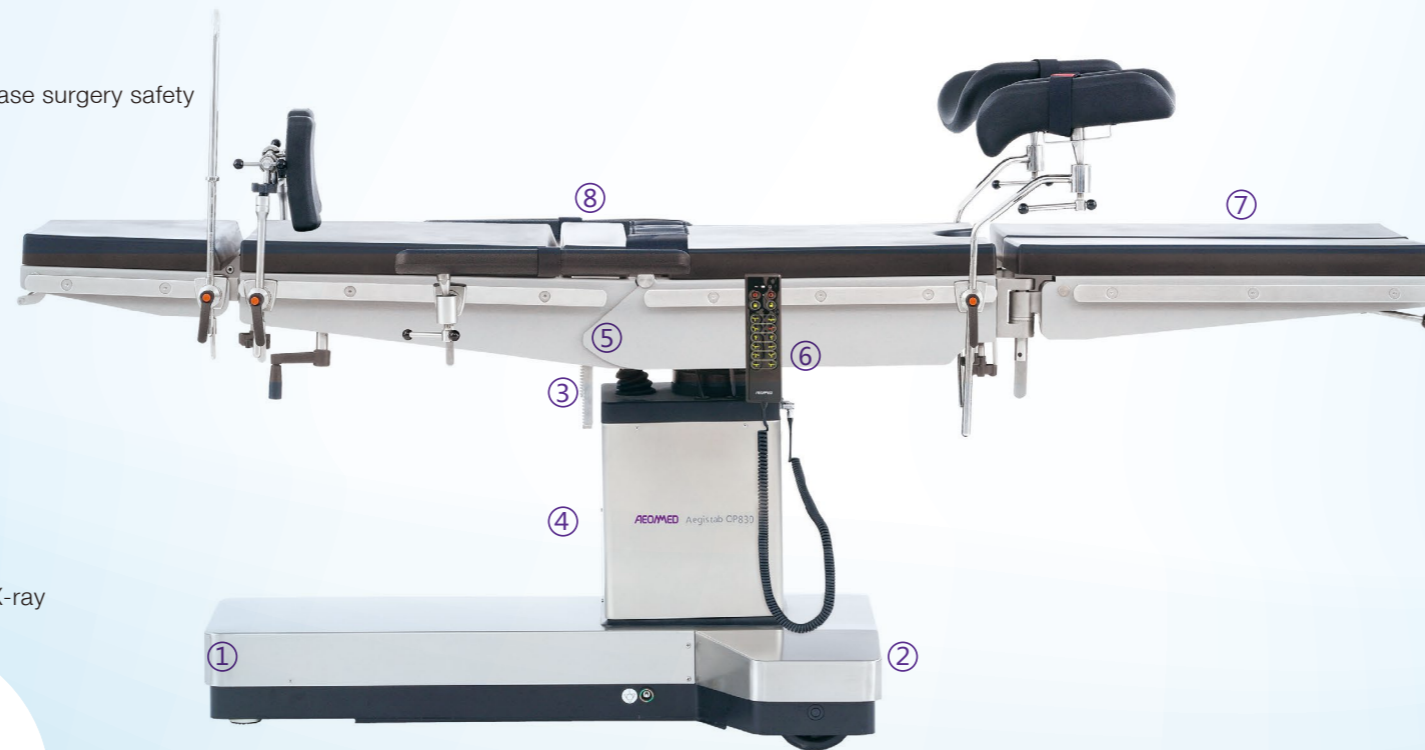


Wide adjustment range for all positions, satisfy demand for different surgeries

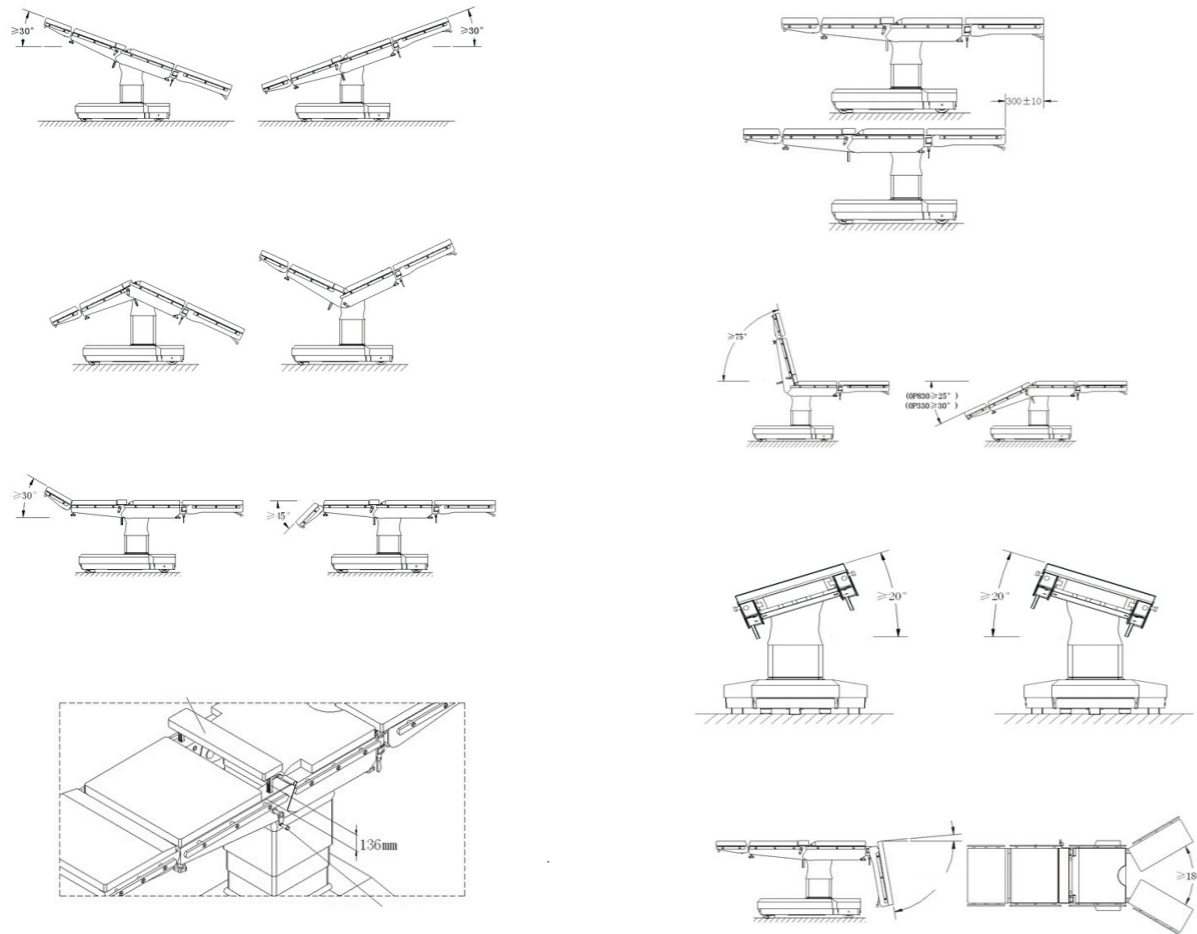

**Option Accessories:
Features:**

All table accessories are made of SUS304 Stainless steel, all surface is easy to clean and disinfect, 6 type standard and 13 option accessories can support a multitude of body positions and procedures, fully meet most of surgeries demands.

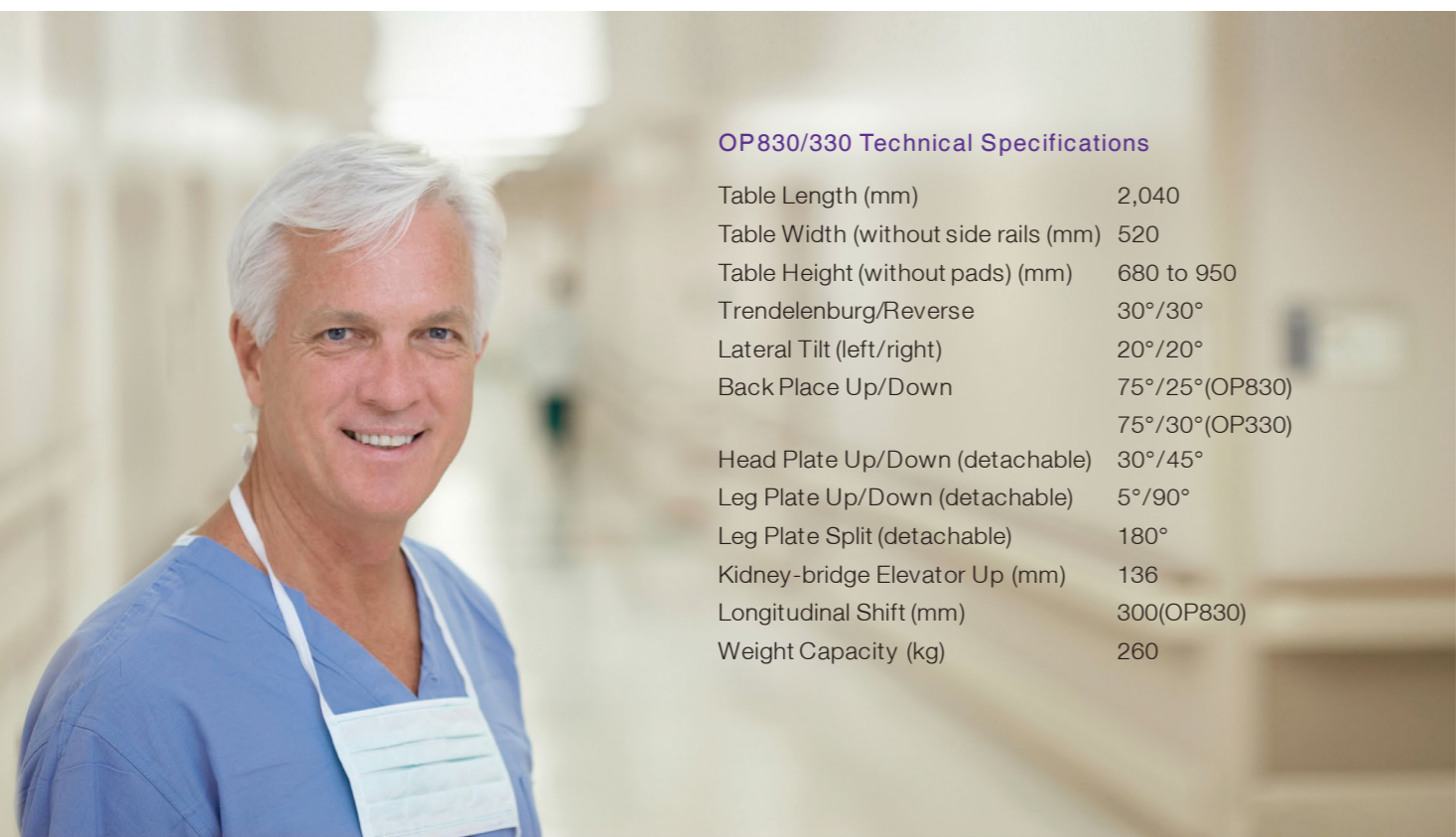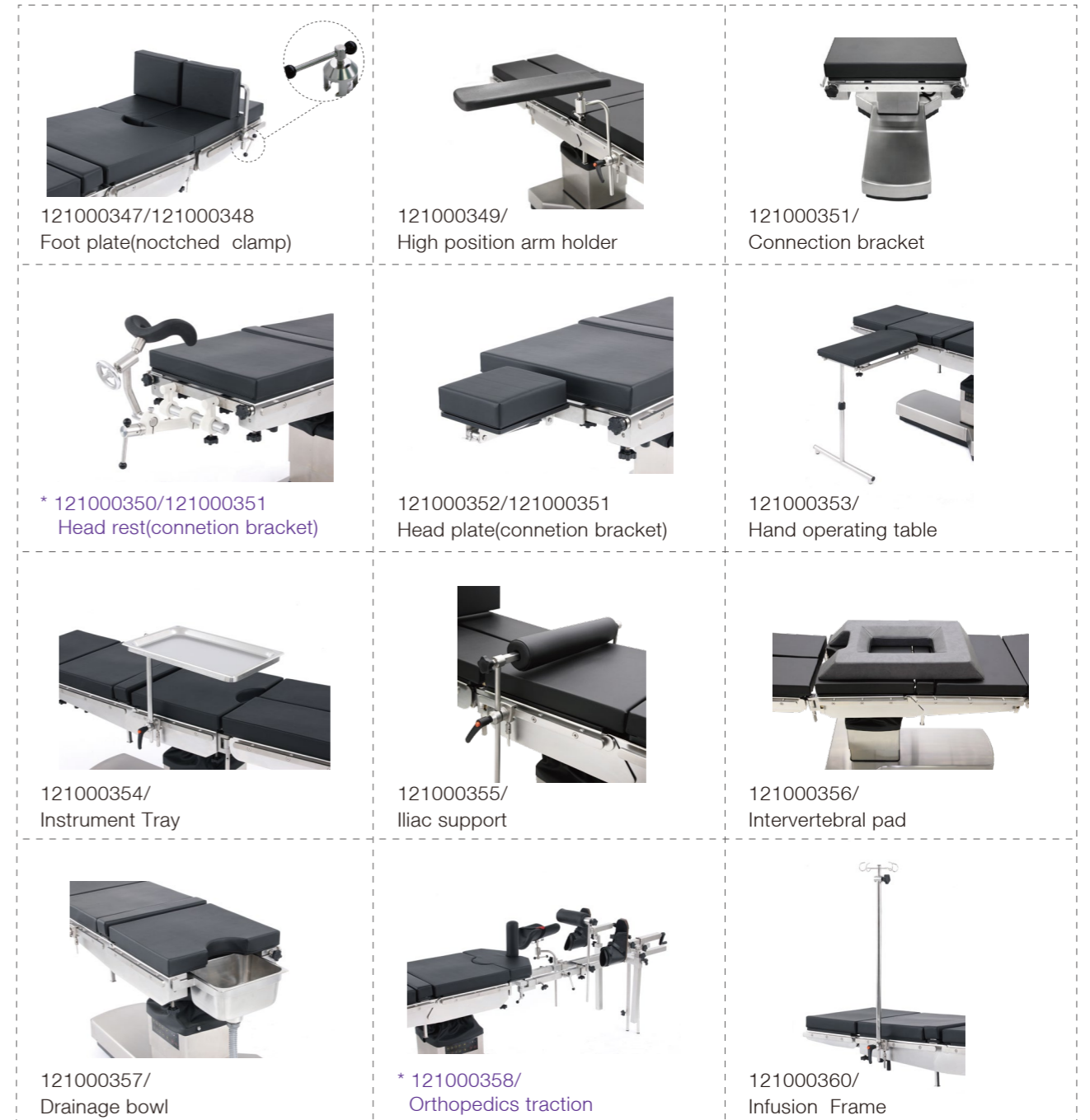


OP830/330 Technical Specifications

Table Length (mm)	2,040
Table Width (without side rails) (mm)	520
Table Height (without pads) (mm)	680 to 950
Trendelenburg/Reverse	30°/30°
Lateral Tilt (left/right)	20°/20°
Back Place Up/Down	75°/25°(OP830) 75°/30°(OP330)
Head Plate Up/Down (detachable)	30°/45°
Leg Plate Up/Down (detachable)	5°/90°
Leg Plate Split (detachable)	180°
Kidney-bridge Elevator Up (mm)	136
Longitudinal Shift (mm)	300(OP830)
Weight Capacity (kg)	260

Orthopedics Features

The Orthopedics Traction is designed based on EU's requirements, meet most of orthopedics surgeries requirement, All part is made of SUS304 stainless steel with high strength, stable during the surgeries. This structure suitable for X-Ray and C-Arm.


Standard Accessories:

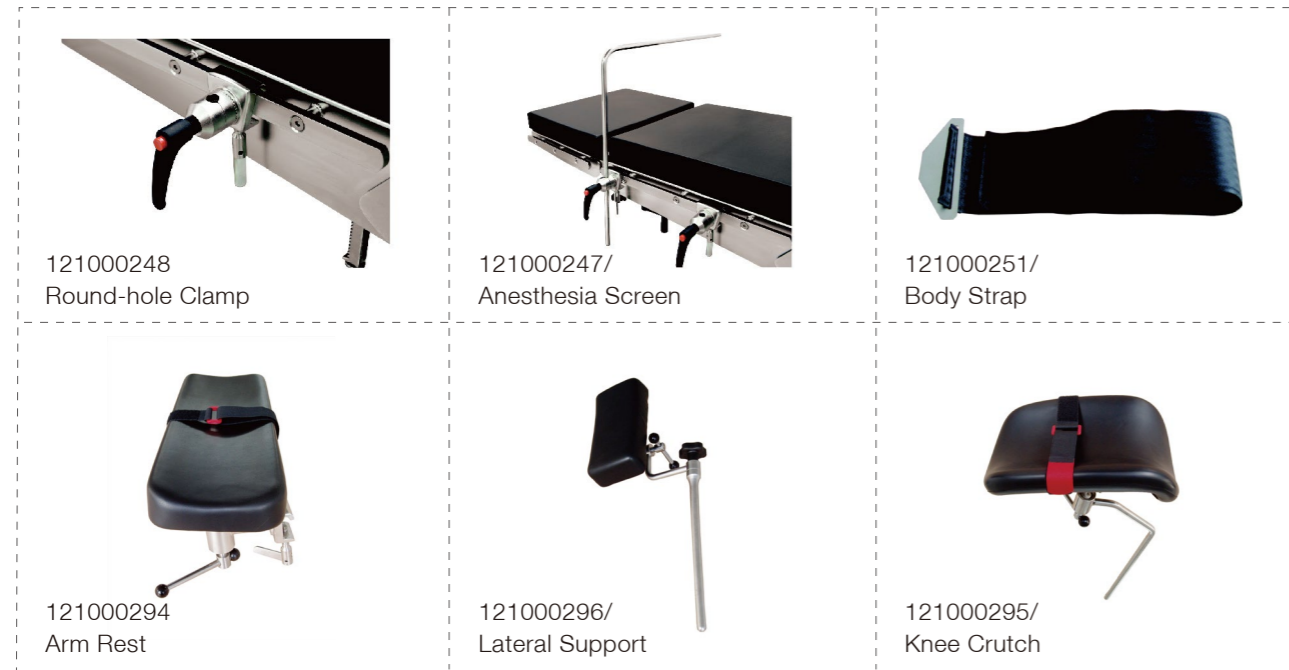

Gel Pad

Functions & Features:

Supplying stable and comfort positioning, Made from Viscoelastic Polymer Gel
Dispersion of pressure, avoiding bedsores Reducing operation time Easy clean out, X-ray through Work temperature from -10°C to 50 °C Good Bio-compatibility with human tissues Cost effective and durable

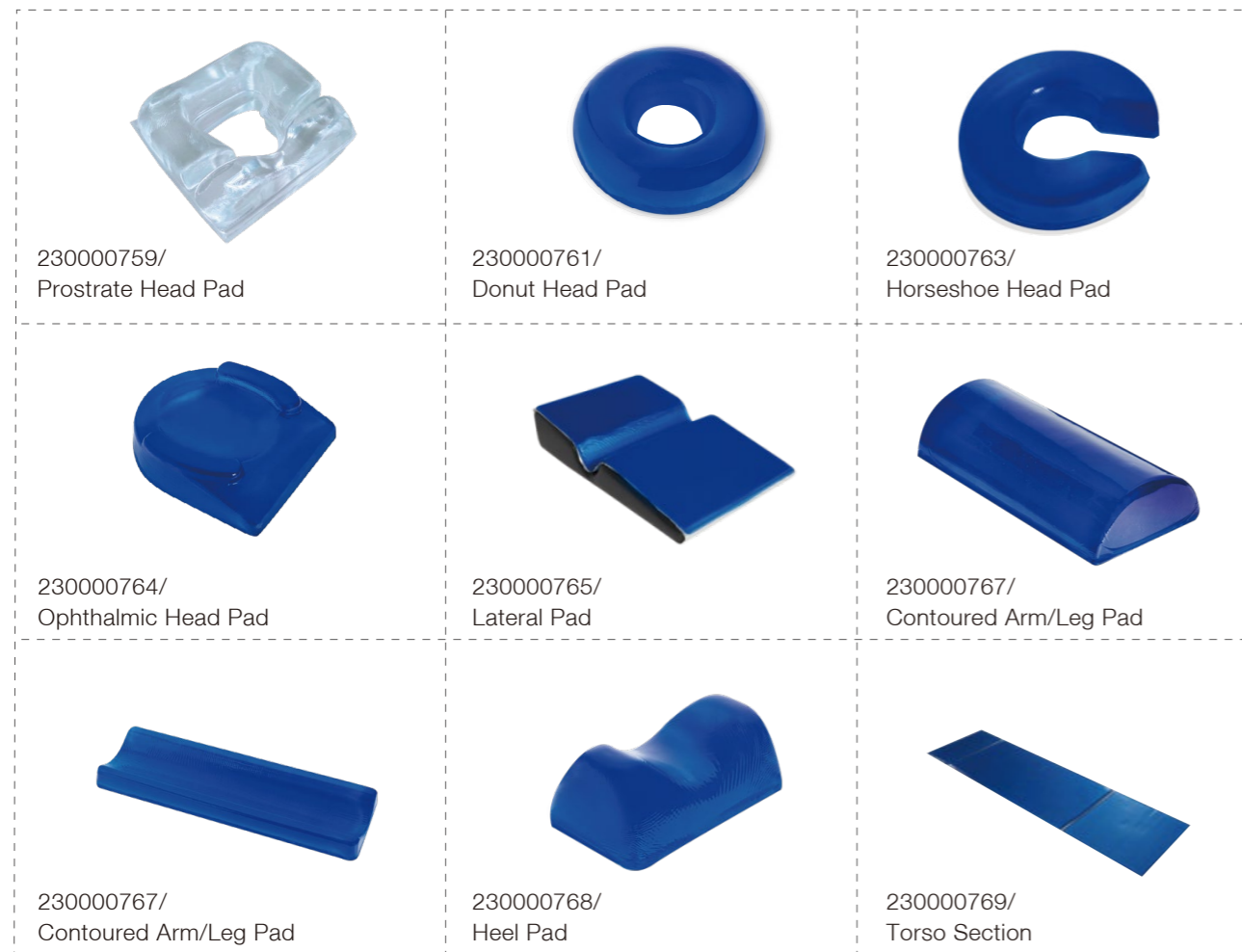

Surgery Positions:



Chest/Abdominal Surgery

Recommended accessories:
230000760 donut head gel pad x1
230000769 Torso Section gel pad x1
121000347 Foot plate 1pair
121000348 Notched Clamp x2
121000294 Arm rest x2



E.N.T Surgery

Recommended accessories:
230000760 donut head gel pad x1



Urology/OB/Gynecology Procedures Surgery

Recommended accessories:
230000760 donut head gel pad x1
230000769 Torso Section gel pad x1
121000295 Knee crutch x2
121000357 Drainage bowl x1
121000294 Arm rest x2



Spinal Surgery

Recommended accessories:
230000759 Prostrate Head Gel Pad x1
121000356 Intervertebral pad x1
121000294 Arm rest x2
230000769 Torso Section gel pad x1
230000766 Universal Dome Shaped Pad x1



Shoulder Operation

Recommended accessories:
121000294 Arm rest x2
121000350 Head rest x1
121000351 Connection bracket x1



Neurosurgery

Recommended Accessories:
121000294 Arm rest x2
121000350 Head rest x1
121000351 Connection bracket x1
230000769 Torso Section gel pad x1
230000766 Universal Dome Shaped Pad x1